


Telefolding and Telescopic Canes



Telefolding Canes Series TF

Telescopic Canes Series TS

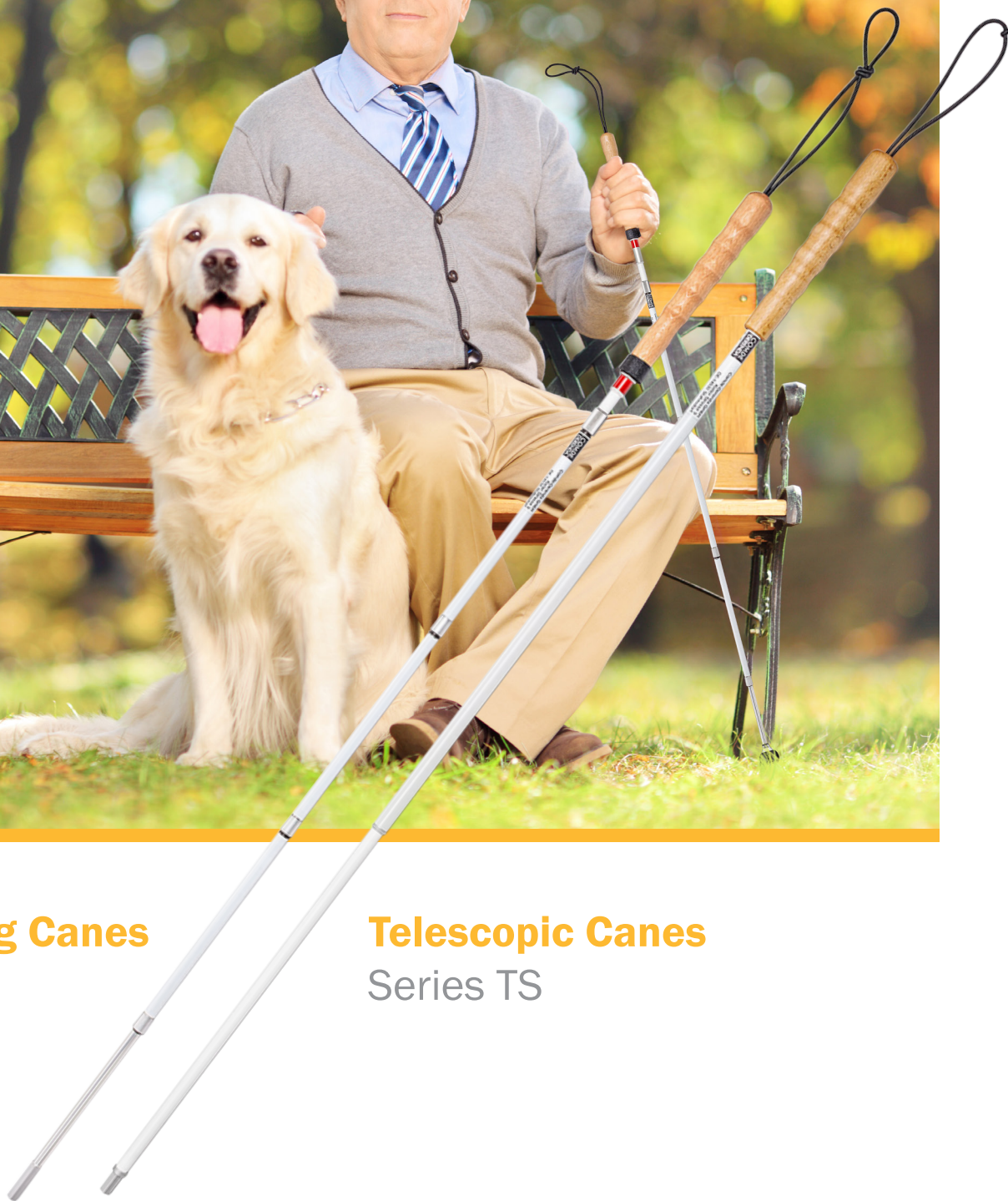

Telefolding Canes

Page 2/4

The Comde-Derenda Telefolding Cane - a groundbreaking design

All Comde-Derenda telefolding canes are made from high-tech aluminum and come with a skin-friendly, ergonomically-shaped beechwood handle.

All canes are length-adjustable and feature a light-reflecting surface. (Prices incl. VAT, excluding cane tips and optional storage pouch)




The multi-segment telefolding cane was developed in cooperation with the Swiss National Association of and for the Blind (SNAB). Thanks to its rugged design, it is practically unbreakable when used normally. The connections between the individual cane segments are free from play due to optimized conical connectors. The connector elements are hard-anodized and therefore very wear-resistant. This design makes the cane easy to use and gives it superior stability and flexural strength.

Thin-walled tube segments guarantee a very long weight. When the cane is extended, its segments are held together by an elastic cord that runs through the inside of the cane.

Before using the cane, the Velcro strip and elastic loop which hold the folded cane together must be unfastened. Then the cane is simply held at the handle as its segments drop down and the cane assembles itself.

The number of cane segments can be selected according to the user's height. Additionally, each cane features a telescopic segment below the handle that enables the user to precisely adjust the cane length. The telescopic segment is loosened and fixed by twisting the grooved tube section below the handle.

The handle features an elastic loop and a Velcro strip. All Comde-Derenda cane tips (or other tips with an M8 thread) can be used with our telefolding canes.

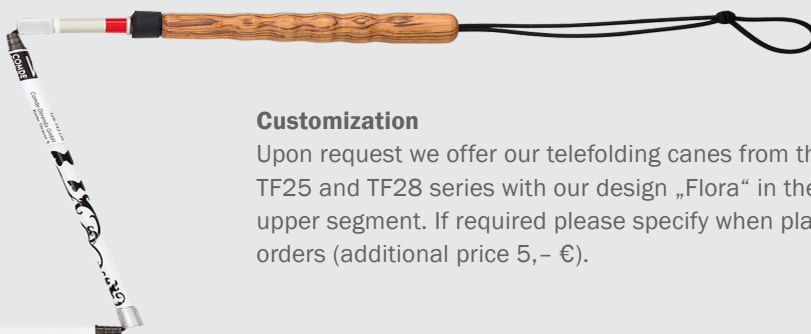

Customization

Upon request we offer our telefolding canes from the TF25 and TF28 series with our design „Flora“ in the upper segment. If required please specify when placing orders (additional price 5,- €).



Telefolding Canes

Page 3/4

Telefolding Canes Series TF25

TF25K



TF25



TF25/28 Hybrid



TF25XL



TF25K

Segments: 5
Length: 88 – 102 cm
Weight: app. 167 g
Item No.: 200004
Price: **174,- €**

TF25

Segments: 6
Length: 111 – 125 cm
Weight: app. 197 g
Item No.: 200003
Price: **189,- €**

TF25/28 HYBRID

Segments: 6
Length: 117 – 134 cm
Weight: app. 205 g
Item No.: 200007
Price: **189,- €**

TF25XL

Segments: 7
Length: 133 – 147 cm
Weight: app. 210 g
Item No.: 200068
Price: **219,- €**

Telefolding Canes Series TF28

TF28K



TF28



TF25/28 Hybrid XL



TF28XL



Smaller sizes available upon request

TF28K

Segments: 5
Length: 100 – 117 cm
Weight: app. 178 g
Item No.: 200002
Price: **174,- €**

TF28

Segments: 6
Length: 126 – 143 cm
Weight: app. 217 g
Item No.: 200001
Price: **189,- €**

TF25/28 HYBRID XL

Segments: 7
Length: 140 – 157 cm
Weight: app. 220 g
Item No.: 200300
Price: **219,- €**

TF28XL

Segments: 7
Length: 150 – 168 cm
Weight: app. 225 g
Item No.: 200000
Price: **219,- €**

Telescopic Canes

Page 4/4

The Comde-Derenda Telescopic Cane - improved reliability

All Comde-Derenda telescopic canes are made from high-tech aluminum and come with a skin-friendly, ergonomically-shaped beechwood handle. All canes are

length-adjustable and feature a light-reflecting surface. (Prices incl. VAT, excluding cane tips and optional storage pouch)

Telescopic Canes Series TS

TS138



TS150



TS138

Segments: 2
Length: 73 – 126 cm
Weight: app. 159 g
Item No.: 200006
Price: 105,- €

TS150

Segments: 2
Length: 86 – 151 cm
Weight: app. 178 g
Item No.: 200005
Price: 105,- €

The telescopic cane is frequently used in mobility training. It is very rugged and will not break during normal use. The cane consists of two anodized aluminum tube segments with telescopic functionality along the entire length of the lower segment. It is very easy to use and its design ensures exceptional stability and flexural strength. Thin-walled tube segments guarantee an extremely low weight.

The upper tube segment features light-reflecting foil, the lower segment is powder-coated. The beechwood handle steamed and waxed for improved comfort. The handle also features an elastic loop.

The cane's length can be precisely adjusted to the user's height by pulling out the lower segment. The cane is then fixed at the desired length by twisting the segment until the built-in expansion plugs lock it in the current position. The telescopic cane is available in two sizes.

All Comde-Derenda cane tips (or other tips with an M8 thread) can be used with our telescopic canes.

Ed. 2019-10